

Whisper

GLOW

English

Model:ND429473



Please read this manual carefully before using your lens and retain it for future reference.

www.thecoolest.com
Printed in Australia

Copyright © 2024 thecoolest Inc.
All Rights Reserved.

The
C  **LEST**



BEFORE WEARING

Wear contact lenses according to the duration prescribed by your prescription, and avoid wearing them for an extended period.



Do not use contact lenses if the sterile packaging has been opened or the protective seal is damaged.



When handling contact lenses, the cleanliness of your hands is of utmost importance. Please maintain good hygiene habits.



Avoid using oily cosmetic products, soaps containing creams, toners, or lotions before wearing your contact lenses.



Check the lense, it should look like a bowl with straight edges. If the edges flair out, it's inside out.




Always remember to follow the advice and instructions of your optometrist regarding the proper handling, wearing, removal, cleaning, disinfection, storage, and use of your contact lenses.




PUTTING ON CONTACT LENSES

Follow along with this quick tutorial to see how to put your lenses in.

With the contact lens on the tip of your index finger, hold your upper and lower eyelids open.

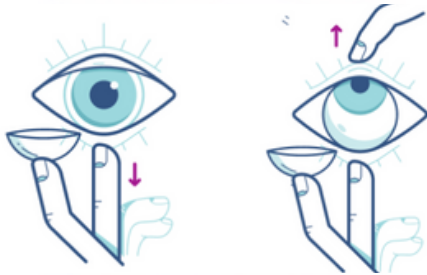

Look up and gently place the lens on your eye. Release the eyelid and close your eye to let the lens settle.

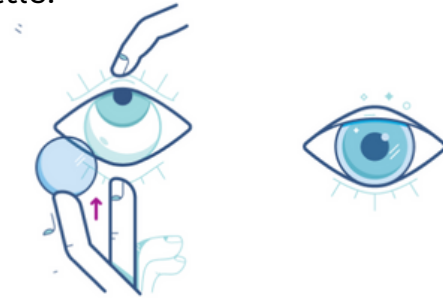

TAKING OFF CONTACT LENSES

Follow along with this quick tutorial to see how to take your lenses out.

Touch the lower edge of the lens with your index finger and slide it down to the lower white part of your eye.



Using your thumb and index finger, lift the lens from your eye.




If you wear bi-weekly or monthly contacts, it's time to clean and disinfect. Daily disposables can be thrown away.



HOW TO CLEAN YOUR CONTACTS

CLEANING CONTACT LENSES

- Remove contact lenses and place them in a clean palm with a few drops of contact lens care solution.
- Gently scrub both sides of the contact lens with your index finger, avoiding excessive force that could damage the lens.
- Carefully rinse the contact lenses, making sure to remove all the care solution and impurities.

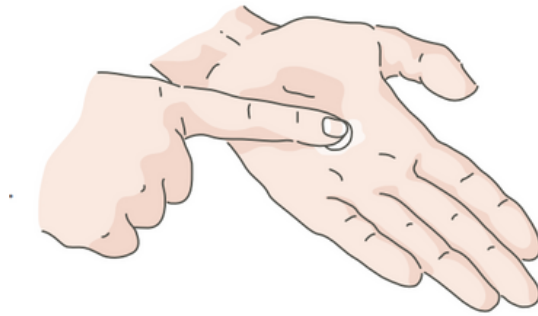

CONTACT LENS DISINFECTION


- Soak your lenses for at least four hours or overnight using a contact lens disinfectant solution to ensure that bacteria and other microorganisms are killed.
- Replace the case at least once a month and rinse the case with care solution before each use.
- Place sterilised contact lenses in a case that has been filled with fresh care solution, making sure that the lenses are completely submerged in the liquid.



HOW TO STORE YOUR CONTACTS

STORING CONTACT LENSES


- To store lenses, disinfect and leave them in the closed/unopened case until ready to wear.
- If you do not wear your lenses immediately following disinfection, you should consult the solution package insert or your Eye Care Practitioner for information on storage of your lenses.
- Always keep your lenses completely immersed in a recommended disinfecting solution when you are not wearing your lenses.

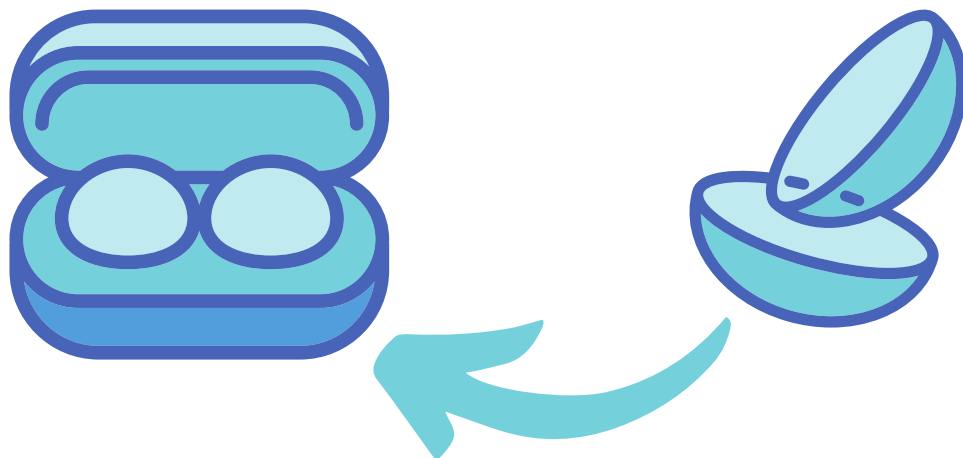

POTENTIAL ADVERSE REACTIONS


When using the Vision Filter Contact Lenses, some users may experience the following adverse reactions:

Visual Illusions: Due to the filter effects, users might experience misalignment of scenes or distortion of object shapes, especially noticeable when first using the lenses. This may lead to errors in judgment during everyday activities such as walking or using transport. It is advised that new users gradually get accustomed to the filter effects in a safe environment.

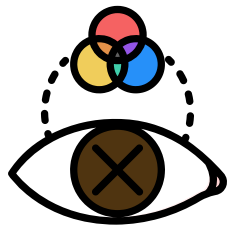

Color Recognition Difficulties: The filters may alter color perception in the environment, making it difficult for users to accurately identify traffic lights or visual cues in emergency situations. Extra attention should be paid to color-related tasks while using these lenses.

A very small number of users may be sensitive to the materials of the lenses or the components of the cleaning solution. If you experience severe allergic symptoms such as significant swelling, pain, or increased discharge from the eyes, stop using the lenses immediately and consult an ophthalmologist.

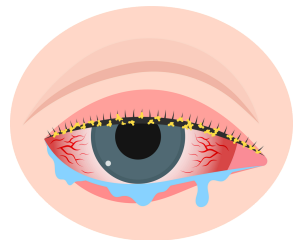

If the above issues occur, **discontinue use** immediately and reassess suitability for continued use of the product. Always follow the guidance in the product manual to ensure safe use.

TROUBLE SHOOTING

Before using contact lenses, users are advised to try them on in non-hazardous, non-emergency situations to get used to the visual changes caused by the filters. When the problems described below are encountered, discontinue use immediately and reassess the suitability of continued use of the product. Always follow the instructions in the product's instruction manual to ensure safe use.

CONTRAINDICATIONS

If you have any the following symptoms, we recommend you NOT TO USE our products

- Severe dry eye syndrome
- Eye infection such as conjunctivitis or keratitis
- Allergic reactions
- Serious eye injury
- Corneal abrasion
- Eye surgery recovery
- Eye inflammation or iritis

These are additional conditions that may impact the suitability of contact lens wear. It's crucial to prioritize eye health and follow your eye doctor's recommendations for safe and comfortable vision correction options.

COMMON PROBLEMS

Question: Tingling in the eye

- Solution: Check the lenses for damage and replace them immediately if they are damaged. And limit the wearing time per day to avoid excessive eye fatigue.

Question: Sensation of a foreign body in the eye

- Solution: Try switching to a different brand or type of lens solution. Alternatively, consult an ophthalmologist to see if there is a need to change the material of the contact lenses.

Question: Dryness in the eyes

- Solution: Prolonged wearing of the lenses may cause dryness. You can use artificial tears or contact lens care products to provide additional moisture.

Question: Recommended wearing time

- Solution: It is recommended not to wear the contact lenses for more than 10 hours. Prolonged wear can cause dryness and blurry vision. If discomfort occurs, remove the contact lenses immediately and seek help from an ophthalmologist.

Question: Lens falling out

- Solution: Ensure that the contact lenses are properly moisturized and correctly worn. Also, try to avoid vigorous activities or rubbing your eyes.

Question: Dizziness or headache after wearing

- Solution: Check if the power of the contact lenses matches your prescription. Ensure that the lenses are properly placed and not inverted. If the issue continues, please consult an ophthalmologist or get in touch with our customer service.

Question: Difficulty removing or the contact lenses shifting

- Solution: Add a few drops of lens solution to moisten the eyes before removing the lenses again. If the lens is out of place, gently massage the closed eyelid to help it return to position.

Question: Difficulty in color recognition

- Solution: If accurate colour recognition is required, it is recommended that contact lenses are removed immediately. And consider using filter options optimised for colour recognition.

Question: Inconsistent or failed filter effect

- Solution: Ensure that the lenses are properly applied and not causing adverse reactions. Check if the lenses are damaged. If the problem persists, quickly remove the contact lenses and contact our customer service for further assistance.

Question: Visual misalignment caused by the filter in specific scenarios

- Solution: In these environments, try to avoid using filters that affect depth perception or cause significant color changes. Choose filters with lower color variation that are more suitable for the scene, or remove the contact lenses if necessary.
- Solution: Under certain filter effects, the enhanced virtual images may conflict with the real world, leading to confusion between recognizing actual objects and virtual effects. It is recommended to gradually adapt and initially use in safe place .

The COOLEST